
RĘCZNY MIKROSKOP WiFi

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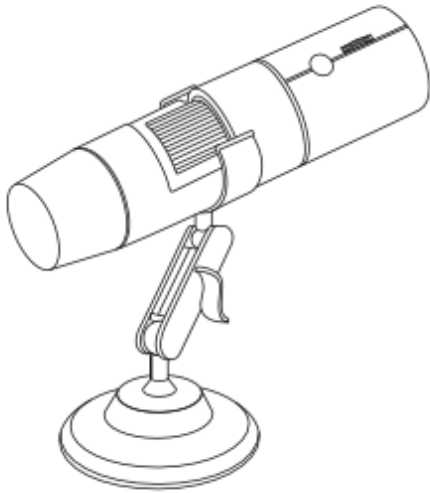
WiFi MANUAL MICROSCOPE

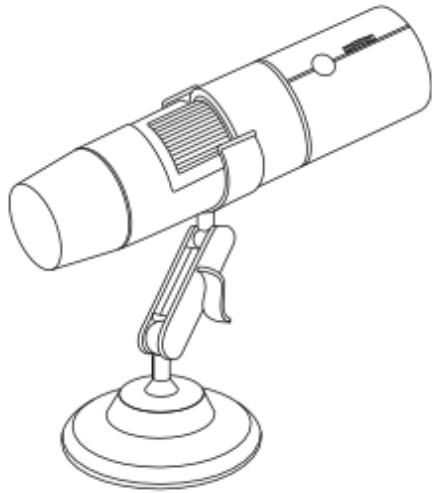Important notes before use


1. Before using the microscope, read the user manual carefully to make sure it will work properly.
2. Fully charge the microscope before using it for the first time. Don't charge directly by PC, you can choose a USB charger (5V 1 A) (on your phone) as it will automatically enter "PC camera mode" when connected to PC. When charging, the LEDs light up red, the charging time is about 2 hours.
3. WiFi only works with the phone, you can't use it on the computer. If you want to use on computer, connect with USB cables and download the corresponding APP software.
4. Before using, remove the lens dust cover on the top of the microscope, otherwise you will not be able to get HD images. Then install the dust cover afterwards to protect the lens from damage and dust.
5. The focus distance of HD image is 3-60mm, adjust the "zoom knob".
6. If you want to observe microscopic worlds such as cells, you need to make a microscope slide yourself and observe it on a mechanical stand.
7. Turn off the useless app on your phone to make sure our microscope runs smoothly.
8. Do not turn off the power when the microscope is connected to a smartphone or computer and do not move it. If there is a sudden cut, the product will be damaged.
9. We are not responsible for any loss or damage to data caused by the use of this product, regardless of the cause or condition of the damage. In addition, we are not responsible for any direct or indirect damage resulting from the failure of this product or its use. Please note.
10. Do not disassemble/repair/recycle/rework this product.
11. Do not use/store this product where there is water, high humidity and a lot of dust and oil smoke.
12. Do not subject this product to strong vibration, impact or excessive force. Also, please do not touch it with wet hands.
13. Do not let children touch this product or accessories unattended.
14. Do not use in unstable places.

15. When cleaning this product, do not use agents containing volatile organic solvents such as benzene and thinner.

16. The product is a simple microscope, the magnification notation and measurement features are only a guideline. It should not be used in cases requiring precise numerical values or in medical facilities related to human life.

17. Do not touch the lens directly with your fingers and do not look directly at the glowing LED.

18. If the digital microscope gets dirty, wipe it with a dry, soft cloth.

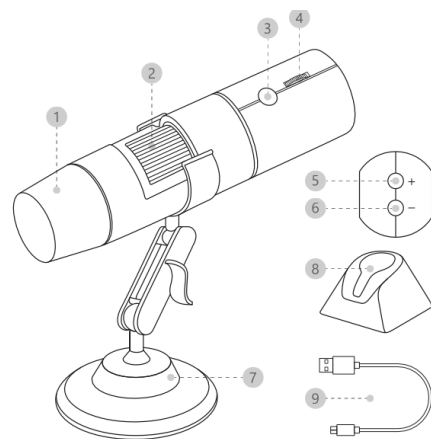Packaging contents:

- 1x WiFi microscope
- 1x Metal microscope holder
- 1x plastic base
- 1x user manual
- 1x USB cable
- 1x user manual
- 1x USB cable

1. Lens protective cover
2. Zoom adjustment wheel
3. Camera power button
4. LED light regulator
5. Zoom in
6. Zoom out
7. Adjustable bracket
8. Base
9. USB cable

Support for Android 6.0 and iOS 9.0 or newer

1. Downloading and installing the app

For iOS: search for the Max-See app in the Apple App Store to download and install it.

For Android: search for the Max-See app in the Google Play store to download and install it.

Note: please follow the installation instructions on your phone and make sure all features are allowed by your phone. Please do not worry that it is security software.


2. Make sure you set it up in this order.

- Download the max-see app, but don't open it.
- Press and hold the power button to turn on the microscope.
- Connect to the max-see app via Wi-Fi in your phone's Wi-Fi settings.
- Open the app.
- Remove the protective cover from the top of the microscope.
- Adjust the LED light using the LED light adjustment knob on the microscope.

Setup instructions

1. First, position the microscope and the object you want to view within sight (3-60 mm).
2. Use the zoom control knob to slowly adjust the multiples. Note: There are 2 best multipliers (50/1000 times) to get clear images.
3. You can also adjust the zoom adjustment button, there are five levels available from multiples of 1.2 to multiples of 2

Application interface instructions

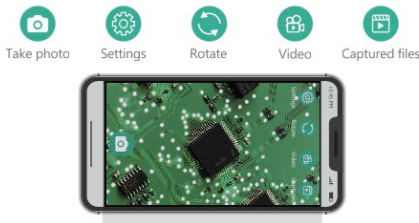
Use on a MAC computer

1. For Macbook (Mac OSX 11.0+).
- Please use the Photo Booth or Quick Time Player software bundled with your Mac.
2. Press and hold to turn on the device and connect it to your MAC with a USB cable.
3. Open Apple's "Photo Booth" application.
4. Select GENERAL-UVC#2 camera options (default), and then start using.

Use on a Windows computer

1. On Windows 10/11, you can connect the device directly to the USB port and search for the phrase CAMERA in WINDOWS, then click it.
- For WIN 7/8 systems, you can solve the problem by downloading Amcap.
- Note: Please disable your laptop's default camera in Windows! And you need to change the privacy settings for the camera, which you need to change to allow access.

2. Or you can download the software from www.inskam.com/download/camera.zip and install Smart Camera.
3. Connect to your computer with a USB cable and long press to turn on the device. The first time you connect the computer, the USB driver will automatically install. Wait for the computer to display a successful installation. If it fails, reconnect it or change the USB port.
4. Click the smart camera icon to start using the microscope.
5. Click the setting options and select the camera in GENERAL-U.

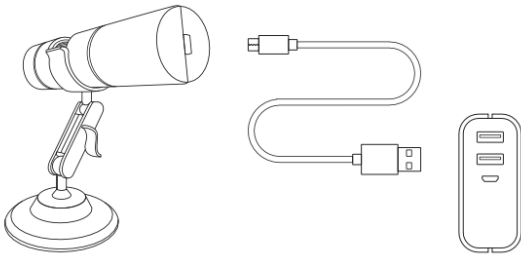Charging

1. Charge the camera fully the first time, do not charge it directly with the computer, choose a 5V 1A USB charger (such as a phone), as it will automatically switch to "PC camera mode" when connected to the computer.
2. When the battery is low, the red indicator light on the back of the device will flash, and the device will automatically turn off after about 10 seconds.
3. While charging, the red indicator light is lit. Charging is complete and the red light is off.
4. The built-in battery is a 900 mAh polymer battery. Normally it takes 3-5 hours to fully charge, it can last for 3 hours.

Product specifications

Matrix resolution: 2MP
 Magnification multipliers: 50-1000X
 Video resolution: 1920*1080P/1280*720P/640*480P
 Video format: MP4/AVI
 Focus distance: 3~60mm (manual adjustment)
 Wi-Fi connection: 10m (open space)
 Image format: JPG
 Computer interface: USB2.0
 Power supply: USB (5V DC)
 Operating system: Android 6.0/IOS9.0/Windows Vista/7/8/10/11
 /MacOS X11.0 or later

Troubleshooting

1. No images:
 - a: Check that the WiFi network is properly connected.
 - b: Close the application and restart.
 - c: Close the microscope and reconnect.
 - d: Check that the Max-See app can be launched at the phone and connect to the Wi-Fi network.
2. LEDs are off, dark or flickering: use the adjustment knob or charge on time.
3. The device malfunctions: press the power button and the zoom button simultaneously for 5 seconds to restore factory reset.
4. Software shuts down automatically: in Apple IOS, find Max See in settings to access photos.

WARNING! THE DEVICE MUST NOT BE DISPOSED OF WITH HOUSEHOLD WASTE.

This symbol indicates that this product must not be disposed of with household waste, according to the EU and your national law. In order to prevent potential damage to the environment or health, the used product must be recycled. In accordance with current legislation, unusable electrical and electronic devices must be collected separately at the designated facilities for recycling, acting on the basis of applicable environmental standards.

Information about used electrical and electronic equipment

The main goal of European and national law regulations is to reduce the amount of waste produced from used electrical and electronic equipment, to ensure an appropriate level of collection, recovery and recycling of used equipment, and to increase public awareness of its harmfulness to the environment, at each stage of use of electrical and electronic equipment. Therefore, it should be pointed out that households play a key role in contributing to reuse and recovery, including recycling of used equipment. The user of electrical and electronic equipment – intended for households – is obliged to return it to authorized collector after its use. However, it should be remembered that products classified as electrical or electronic equipment should be disposed of at authorized collection points.